

Beyond i700's Proven Tech

Still the same **super fast**, **light**, and **accurate** *i***700**. But now **wireless**.



No More Wires

Free from wires, scanning is comfortable at any angle. Scan proximal areas effortlessly without constraints in movement.



Fast, Accurate. Connected

Scanning speed of up to 70 FPS*. You can experience the amazingly smooth and fast scanning performance of Medit's best-selling model. Enjoy the fast speed of the *i*700 wirelessly, flawlessly.



Supercharged Yet Light


Outstanding weight balance even with a battery attached for added comfort.



Powerful Hardware Meets Versatile Software

Optimized for convenience. Not only from hardware aspects, but also software. *i*700 wireless, packed with the latest of MEDIT technology.

*Frames per second: number of frames of scanning data per second

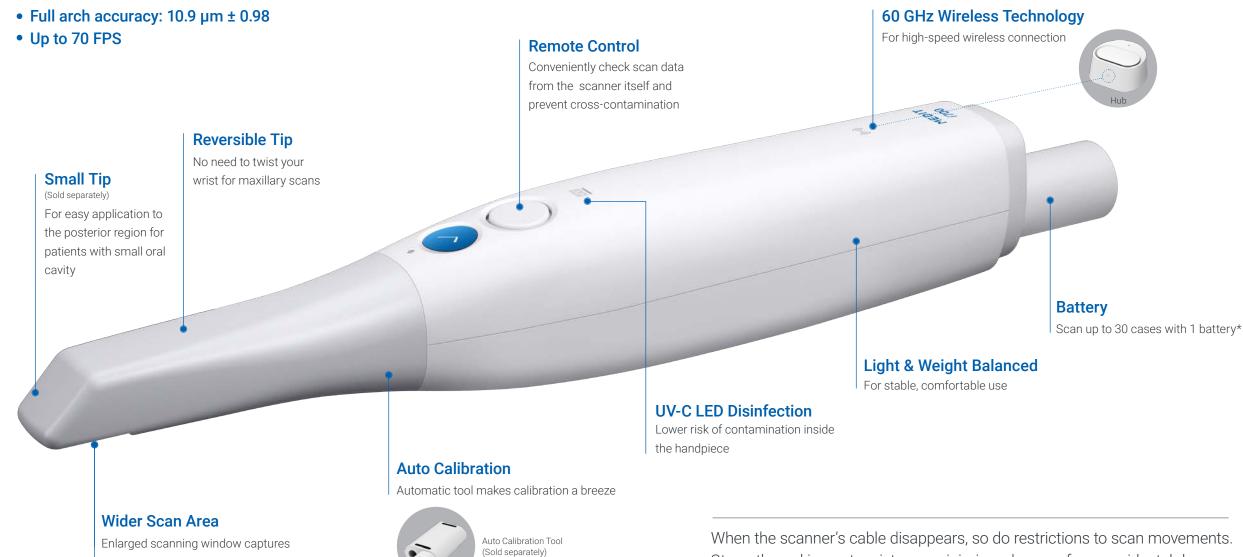

large areas at once while maintaining

the small size of the scanner head


Strengthened impact resistance minimizes damage from accidental drops. Weight balanced design of the *i*700 wireless provides a stable, comfortable scanning experience.

*Performance may vary according to usage conditions.



The Best Battery Life Ever

1 hour of continuous scanning, up to 8 hours in standby mode

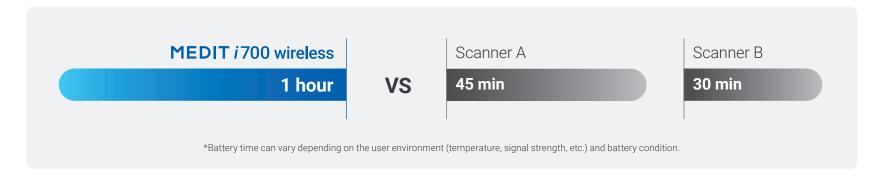

All-Day-Scan-with-Single-Battery. Scan all day without worrying about charging.




Maxillary + Mandibular + Occlusion

Total scan time of

2 mins*



Scan
up to **30 cases**with 1 battery

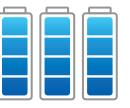

Up to **90 cases** using 3 included batteries



Intelligent power management function automatically switches the device to sleep mode when not in use, standing by up to 8 hours.

It takes just seconds to wake the device for scanning.



Wired power connection:

Once you've used up all 3 included batteries, simply connect the power cable and keep using the device.

Magical Scanning, **Easier Than Ever**



Experience amazing digital dentistry with **Medit Scan for Clinics**. Our cutting-edge software technology makes scanning more convenient for you.


For scanning assistance

Smart Scan Guide



Strengthened user guidance provides Enhanced A.I. filtering performance help for abnormalities detected during eases the scanning process by filtering such as the color of medical gloves, scanning.

Smart Scan Filtering

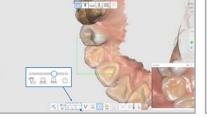


out soft tissue such as lips, cheeks, and tongue.

Smart Color Filtering



Automatically filters registered colors,

Smart Stitching

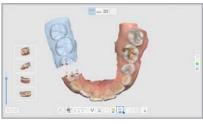

Scan with your own scanning strategy Fragmented scan data are aligned with automatically detected common area.

For more precise and accurate scan data

Smart Shade Guide




Automatically distinguishes and suggests the best matching shade from scan data.


A.I. Scan Body Matching



Secures accurate scan body data in Secures accurate margin data in difficult-to-scan areas, such as the difficult-to-scan situations posterior region.

A.I. Abutment Matching

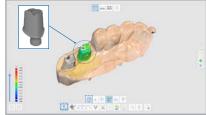

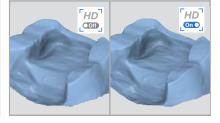

Get a sharper scan image in the margin area with the upgraded HD Mode.

Medit Link

Begin your clinic's data management with Medit Link.



You can create and manage scan data matched with patient information in the electronic chart program and link them with treatment solutions of various workflows. The various apps of the Medit Link App Box help provide higher standard consultation that is easier to understand for patients.

Medit Link's cloud storage encrypts and stores data for reliable protection against virus and ransomware attacks. It also enables data sharing and real-time communication between dentists, patients, and lab technicians anytime, anywhere on a PC, a tablet, or any mobile device with an internet connection.

Digital Dentistry with **Medit Link** & **Medit Apps**

Digitize Your Treatment with Medit Link



Use Your Scan Data for Much More

Scan data created with the Medit *i*700 wireless can be linked to the workflow of various solution providers via Medit Link. The data can also be used with any software that can open standard 3D formats such as STL, PLY, and OBJ. And the wide utilization of our open system is expanding as we speak.

What you can do in Medit Link:


Clear aligners | Crown and bridges | Implant-supported restorations | Full and partial dentures | Orthodontics | Implant surgical guides | Sleep apnea devices | Dental models | 3D printing...



Medit Apps

Medit Apps quickly introduce digital workflows using scan data and help you to utilize them in a diverse range of cases. You can analyze data, create restorations and models, and communicate with your patients in the simplest way. With Medit Apps, a lot more can be done on the chairside like designing sprints, temporaries, bridges, and so on! By the way, you don't need to have a Medit scanner or pay a huge annual fee to join the Medit Ecosystem. Feel literally free to enjoy digital dentistry with Medit Apps.

→ See Next Page for **Medit Ecosystem**

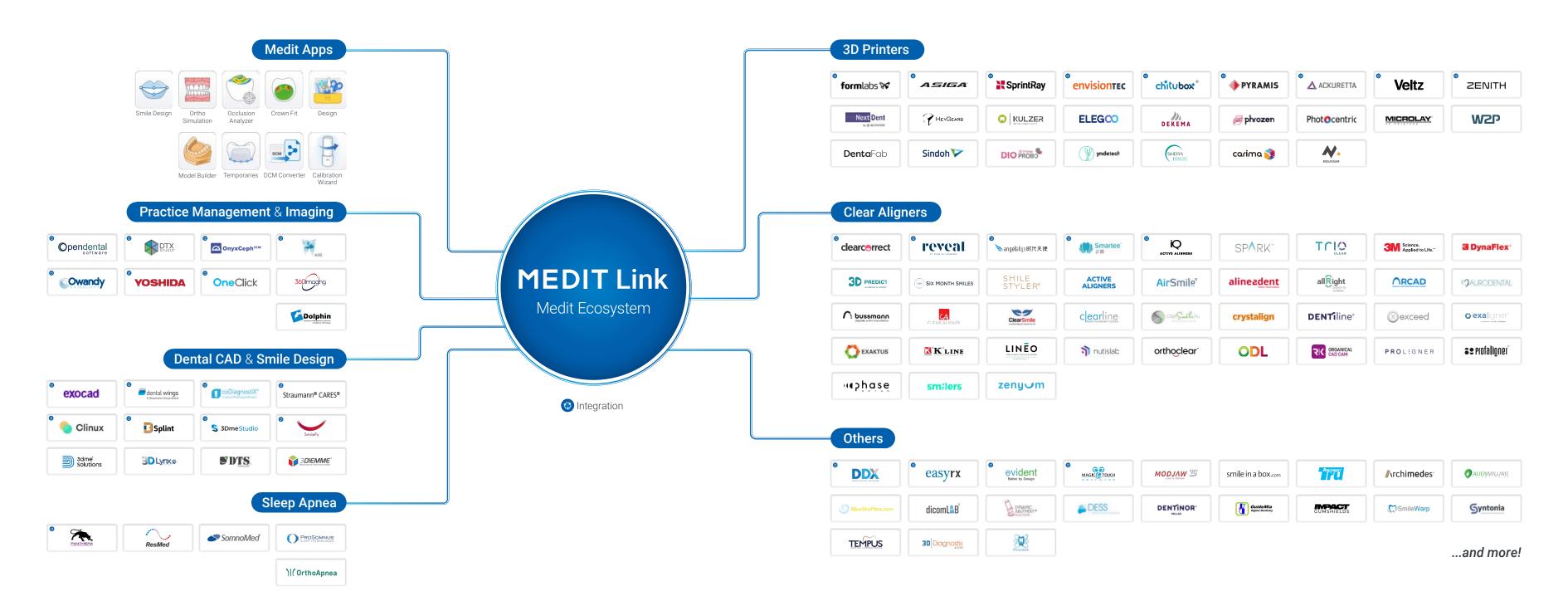

Specifications

Category		Medit <i>i</i> 700 wireless
Scanning Technology	Scanning Frame	Up to 70 FPS
	Imaging Technology	3D-in-motion video technology 3D full color streaming capture
	LED Source	LED
	Anti-Fogging Technology	Adaptive anti-fogging
Accuracy	Full-Arch	$10.9 \mu m \pm 0.98$
Handpiece	Dimensions	313 x 44 x 47.4 mm
	Total Weight	328 g (Includes battery, tip)
Tip	Tip Size	22.5 x 17.1 mm
	Scan Area	15 x 13 mm
	Autoclavable	Autoclave 150 times 121°C 30 mins (Gravity Type) 134°C 4 mins (Pre-Vacuum Type) 135°C 10 mins (Gravity Type)
	Reversible Tip	Yes
Special Mode	Remote Control Mode	Yes
	UV-C Disinfection	Yes
Connectivity		Wireless USB 3.1 Gen1 (Wireless Hub)



About Us

Medit specializes in 3D measurement and CAD/CAM solutions for dental clinics and labs, including intraoral scanners, based on its own patented state-of-the-art technology. We also develop platform solutions for digital dentistry, supporting collaborative workflows. Additionally, the Solutionix line provides 3D scanners and software to the industrial market. Our goal is to provide innovative technology and the highest quality products to ensure mutual growth for all partners.

Our Medit Community







User Group N

Help Cent

Continuous feedback, suggestions of new treatment workflows, and active discussion.

From hot debates to fun chats, the Medit community is always lively with our passionate customers.

Colleagues from all over the world actively support and help one another, and sincerely welcome new

colleagues into digital dentistry.

Scan the QR code to join!

Medit and its community of colleagues support and wish you the best of luck on your new journey into digital dentistry.





English